


Summary

Name	Semenogelin-1/SEMG1
Purity	Greater than 95% as determined by reducing SDS-PAGE
Endotoxin level	<1 EU/µg as determined by LAL test.
Construction	Recombinant Human Semenogelin-1 is produced by our Mammalian expression system and the target gene encoding Gln24-Thr402 is expressed with a 6His tag at the C-terminus.
Accession #	AAH07096.1
Host	Human Cells
Species	Human
Predicted Molecular Mass	43.8 KDa
Formulation	Lyophilized from a 0.2 μ m filtered solution of 20mM Hac-NaAc, 150mM NaCl, pH 4.5.
Shipping	The product is shipped at ambient temperature. Upon receipt, store it immediately at the temperature listed below.
Stability&Storage	Store at \leq -70°C, stable for 6 months after receipt. Store at \leq -70°C, stable for 3 months under sterile conditions after opening. Please minimize freeze-thaw cycles.
Reconstitution	Always centrifuge tubes before opening. Do not mix by vortex or pipetting. It is not recommended to reconstitute to a concentration less than $100\mu g/ml$. Dissolve the lyophilized protein in distilled water. Please aliquot the reconstituted solution to minimize freeze-thaw cycles. Always centrifuge tubes before opening. Do not mix by vortex or pipetting. It is not recommended to reconstitute to a concentration less than $100\mu g/ml$. Dissolve the lyophilized protein in distilled water. Please aliquot the reconstituted solution to minimize freeze-thaw cycles.

SDS-PAGE image

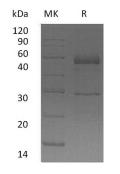

Background



Alternative Names	Semenogelin-1; Semenogelin I; SGI; SEMG1; SEMG; Alpha-Inhibin-92; Alpha- Inhibin-31; Seminal Basic Protein
Background	Semenogelin-1 (SEMG1) is the predominant protein in semen; it is a secretory protein involved in the formation of a gel matrix entrapping the accessory gland secretions and ejaculated spermatozoa. The prostate-specific antigen (PSA) protease processes SEMG1 into smaller peptides, each possibly having a separate function. In the proteolysis process, Alpha-inhibin-92 and alpha-inhibin-31 are produced; they inhibit the secretion of pituitary follicle-stimulating hormone. At the same time, it breaks down the gel matrix, allowing the spermatozoa to move more freely.

Note

For Research Use Only, Not for Diagnostic Use.