

# **Product Name:** cAMP Protein Kinase Catalytic Subunit Rabbit Monoclonal Antibody Catalog #: AMRe04058

For research use only.

### **Summary**

**Description** Recombinant rabbit monoclonal antibody

**Host** Rabbit

Application WB,IHC,ICC/IF,IP

Reactivity Human,Mouse,Rat

Conjugation Unconjugated

Modification Unmodified

**Isotype** IgG

Clonality Monoclonal
Form Liquid

Concentration 0.51mg/ml. The concentration of this product may be batch-dependent.

Storage Aliquot and store at -20°C (valid for 12 months). Avoid freeze/thaw cycles.

**Shipping** Ice bags

50mM Tris-Glycine(pH 7.4), 0.15M NaCl, 40% Glycerol, 0.01% Sodium azide and 0.05% **Buffer** 

protective protein

**Purification** Affinity Purification

#### **Application**

**Dilution Ratio** WB 1:500-1:1000,IHC 1:50-1:100,ICC/IF 1:50-1:200,IP 1:20-1:50

Molecular Weight Calculated MW: 41 kDa; Observed MW: 41 kDa

## **Antigen Information**

Gene Name PRKACA

Alternative Names PKACA; PPNAD4

 Gene ID
 5566

 SwissProt ID
 P17612

**Immunogen** A synthetic peptide corresponding to target protein

## **Background**

Phosphorylates a large number of substrates in the cytoplasm and the nucleus. Regulates the abundance of compartmentalized pools of its regulatory subunits through phosphorylation of PJA2 which binds and ubiquitinates these

Web: https://www.enkilife.com E-mail: order@enkilife.com techsupport@enkilife.com Tel: 0086-27-87002838




subunits, leading to their subsequent proteolysis. Phosphorylates CDC25B, ABL1, NFKB1, CLDN3, PSMC5/RPT6, PJA2, RYR2, RORA and VASP. RORA is activated by phosphorylation. Required for glucose-mediated adipogenic differentiation increase and osteogenic differentiation inhibition from osteoblasts. Involved in the regulation of platelets in response to thrombin and collagen; maintains circulating platelets in a resting state by phosphorylating proteins in numerous platelet inhibitory pathways when in complex with NF-kappa-B (NFKB1 and NFKB2) and I-kappa-B-alpha (NFKBIA), but thrombin and collagen disrupt these complexes and free active PRKACA stimulates platelets and leads to platelet aggregation by phosphorylating VASP. Prevents the antiproliferative and anti-invasive effects of alpha-difluoromethylornithine in breast cancer cells when activated. RYR2 channel activity is potentiated by phosphorylation in presence of luminal Ca2+, leading to reduced amplitude and increased frequency of store overload-induced Ca2+ release (SOICR) characterized by an increased rate of Ca2+ release and propagation velocity of spontaneous Ca2+ waves, despite reduced wave amplitude and resting cytosolic Ca2+. PSMC5/RPT6 activation by phosphorylation stimulates proteasome. Negatively regulates tight junctions (TJs) in ovarian cancer cells via CLDN3 phosphorylation. NFKB1 phosphorylation promotes NF-kappa-B p50-p50 DNA binding. Involved in embryonic development by down-regulating the Hedgehog (Hh) signaling pathway that determines embryo pattern formation and morphogenesis. Prevents meiosis resumption in prophase-arrested oocytes via CDC25B inactivation by phosphorylation. May also regulate rapid eye movement (REM) sleep in the pedunculopontine tegmental (PPT). Phosphorylates APOBEC3G and AICDA. Isoform 2 phosphorylates and activates ABL1 in sperm flagellum to promote spermatozoa capacitation. Phosphorylates HSF1; this phosphorylation promotes HSF1 nuclear localization and transcriptional activity upon heat shock (PubMed:21085490).

#### **Research Area**

**Signal Transduction** 


#### **Image Data**

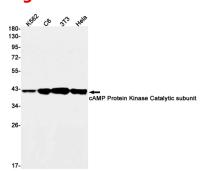


Western blot analysis of cAMP Protein Kinase Catalytic subunit in K562, C6, 3T3, Hela lysates using cAMP Protein Kinase Catalytic Subunit antibody.

Web: https://www.enkilife.com E-mail: order@enkilife.com techsupport@enkilife.com Tel: 0086-27-87002838



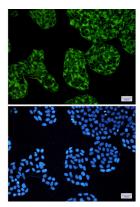

Immunocytochemistry analysis of cAMP Protein Kinase Catalytic subunit (green) in Hela using cAMP Protein Kinase Catalytic subunit antibody,and DAPI(blue)

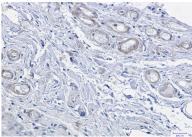

Immunohistochemistry analysis of paraffin-embedded Human colon cancer using cAMP Protein Kinase Catalytic subunit antibody. High-pressure and temperature Sodium Citrate pH 6.0 was used for antigen retrieval.

Web: https://www.enkilife.com E-mail: order@enkilife.com techsupport@enkilife.com Tel: 0086-27-87002838