

Product Name: Kappa light chain (15R1) Rabbit Monoclonal Antibody Catalog #: AMRe12897

For research use only.

Summary

Description Recombinant rabbit monoclonal antibody

Host Rabbit

Application WB,IHC,IP,IF-P

Reactivity Human

ConjugationUnconjugatedModificationUnmodified

Isotype IgG

Clonality Monoclonal

Form Liquid

Concentration 0.5mg/ml. The concentration of this product may be batch-dependent.

Storage Aliquot and store at -20°C (valid for 12 months). Avoid freeze/thaw cycles.

Shipping Ice bags

Rabbit IgG in phosphate buffered saline, pH 7.4, 150mM NaCl, 0.02% New type preservative

Buffer N and 50% glycerol. Store at +4°C short term. Store at -20°C long term. Avoid freeze / thaw

cycle.

Purification Affinity purification

Application

Dilution Ratio WB 1:2000-1:20000,IHC 1:1000-1:5000,IP 1:10-1:100,IF-P 1:1000-1:5000

Molecular Weight 12kDa

Antigen Information

Gene Name IGKC

Alternative Names HCAK1; Ig kappa chain C region; IGKCD; Immunoglobulin InV;

Gene ID

SwissProt ID P01834

Immunogen Full length native protein

Background

The five types of immunoglobulin heavy chains are known as: IgG, IgA, IgM, IgD, and IgE. IgG is divided into four subclasses,

Web: https://www.enkilife.com E-mail: order@enkilife.com techsupport@enkilife.com Tel: 0086-27-87002838




and IgA is divided into two subclasses. In serum IgA and IgG are monomers with a single 4 polypeptide unit; while, IgM is a pen tamer. IgA may also form polymers. Kappa light chain antibody can be used for the identification of leukemias, plasmacytomas and certain non Hodgkin's lymphomas. Constant region of immunoglobulin light chains. Immunoglobulins, also known as antibodies, are membrane-bound or secreted glycoproteins produced by B lymphocytes. In the recognition phase of humoral immunity, the membrane-bound immunoglobulins serve as receptors which, upon binding of a specific antigen, trigger the clonal expansion and differentiation of B lymphocytes into immunoglobulins- secreting plasma cells. Secreted immunoglobulins mediate the effector phase of humoral immunity, which results in the elimination of bound antigens (PubMed:22158414, PubMed:20176268). The antigen binding site is formed by the variable domain of one heavy chain, together with that of its associated light chain. Thus, each immunoglobulin has two antigen binding sites with remarkable affinity for a particular antigen. The variable domains are assembled by a process called V-(D)-J rearrangement and can then be subjected to somatic hypermutations which, after exposure to antigen and selection, allow affinity maturation for a particular antigen (PubMed:17576170, PubMed:20176268).

Research Area

Immunology

Image Data

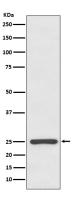

Western blot analysis of Kappa light chain expression in human plasma lysate.

Web: https://www.enkilife.com E-mail: order@enkilife.com techsupport@enkilife.com Tel: 0086-27-87002838